



**K-FILES**



**K-REAMERS**



**H-FILES**

Stainless steel

### Applications:

- K-files: for preparing canals (marking and boring).
- K-reamers: for preparing canals.
- H-files: for scraping walls and removing the resulting debris.

## Why choose R&S K-files, K-reamers and H-files?

- + Easy to identify**  
ISO colour coded.
- + Silicone mounted stops**  
Perfectly adjusted to the canal length.
- + Ergonomic handles**  
Make the tools easier to grip.
- + Safe to use**  
Tip rubber ensures smooth handling.
- + Full range**  
A wide selection of sizes and lengths (15 diameters and 2 lengths).



## Advice from our dental experts

- Rinse thoroughly with a stabilised sodium hypochlorite solution.
- Work in the canal using back-and-forth movements without exerting pressure.
- Do not use deformed instruments in root canal procedures.

### Technical specifications and instructions for use



#### K FILES

A ground blade with a cross-section that is first square then triangular and with a tighter thread than a reamer for both cutting power and flexibility.



#### K-REAMERS

A ground blade with a square cross-section designed for greater cutting power and flexibility.



#### H-FILES

Instrument shaped in the shank by grinding. The thread is progressive to obtain regular flexibility of the instrument on its active part.

#### Description

#### Standard

#### Taper (%)

#### Cross-section shape

#### Materials

#### ISO (Ø)

#### Length (mm)

#### Packaging

#### Sterilisation procedure

#### Number of uses (under normal conditions)

ISO 3630-1



2%

Square for numbers 06, 08 and 10.  
Triangular for numbers 15 to 80.

Square.

Round with just one very sharp cutting edge.

Handle: plastic  
Stop: silicone  
Blade: stainless steel

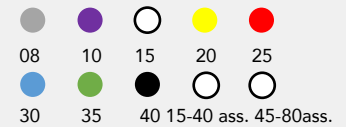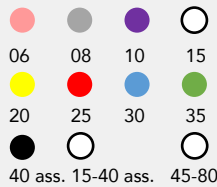

21 or 25

Box of 6 instruments.

Sterilisation in an autoclave at 134°C, for a minimum of 18 min. in sterilisation pouch.

2 for numbers 06 to 20  
5 for the other numbers.

2 for numbers 06, 08 and 10  
5 for the other numbers.

2 for numbers 06 to 20  
5 for the other numbers.

## Additional products

#### Paper points Gutta percha points



#### Nerve broaches



#### Root canal pluggers

