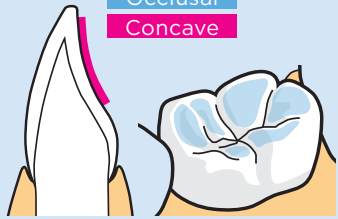
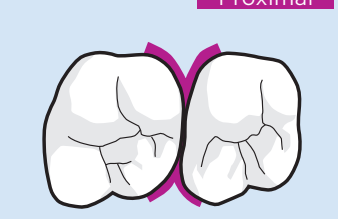
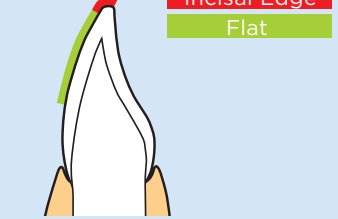
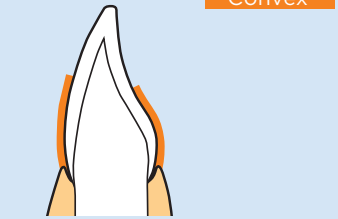





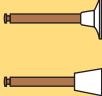



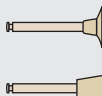
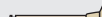
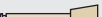


# Enhance® / Enhance® PoGo® Finishing and Polishing Tips



<p><b>Site</b></p> <p><b>Steps</b></p>	<p>Occlusal Concave</p> 	<p>Proximal</p> 	<p>Incisal Edge Flat</p> 	<p>Convex</p> 
<p><b>Gross Reduction</b></p>	<p>Diamond or Carbide bur</p> 	<p>Diamond or Carbide bur</p> 	<p>Diamond or Carbide bur</p> 	<p>Diamond or Carbide bur</p> 
<p><b>Excess Removal + Contouring + Finishing</b></p>	<p>Finishing Point</p>  <p><b>Enhance®</b></p>	<p>Finishing Disc or Finishing Cup</p>  <p><b>Enhance®</b></p>	<p>Finishing Disc</p>  <p><b>Enhance®</b></p>	<p>Finishing Cup</p>  <p><b>Enhance®</b></p>
<p><b>Polishing High Gloss</b></p>	<p>Polishing Point</p>  <p><b>Enhance® PoGo®</b></p>	<p>Polishing Disc or Polishing Cup</p>  <p><b>Enhance® PoGo®</b></p>	<p>Polishing Disc</p>  <p><b>Enhance® PoGo®</b></p>	<p>Polishing Cup</p>  <p><b>Enhance® PoGo®</b></p>

## Impressive Results

## Ordering Information



Class II DO restoration with ceram.x® universal after finishing with Enhance®



Class II DO restoration with ceram.x® universal, final result after polishing with Enhance® PoGo®

### Enhance® Finishing System Refills

40 Finishing Discs  
40 Finishing Cups  
40 Finishing Points

624045  
624055  
624065



### Enhance® PoGo® Polishing System Refills

40 Polishing Discs  
40 Polishing Cups  
40 Polishing Points

662010  
662023  
662022



### Enhance® /Enhance® PoGo® Complete System

Enhance® Finisher  
10 each: Discs, Cups, Points  
Enhance® PoGo® Polisher  
10 each: Discs, Cups, Points

662020



For more information visit [www.dentsplysirona.com](http://www.dentsplysirona.com).

Find us on  [www.facebook.com/dentsplysirona.restorative](https://www.facebook.com/dentsplysirona.restorative)

K79200251-00  
© Dentsply Sirona 5/2017