

3M Science.
Applied to Life.™

3M™ RelyX™ Universal
Resin Cement

Goodbye
complexity.



The two that can do it all.

3M™ RelyX™ Universal Resin Cement works both as a standalone, self-adhesive cement and as an adhesive cement when combined with 3M™ Scotchbond™ Universal Plus Adhesive.

Together they eliminate the hassle and confusion of multiple resin cements, primers and adhesives, **simplifying your direct and indirect restorative workflows.** With fewer products to stock and replenish, you save space, time and money.

“Excellent, hygienic dispensing system with very little waste.”

Dr M.W. (Manchester, UK)***

80%
less cement
waste* with
Micro Mixing Tip

More hygienic.**
Self-sealing for
storage without
used mixing tip.

“Significantly less waste and therefore a much more cost effective system.”

Dr K.S. (London, UK)***

More ergonomic
shape with
50% less
plastic waste*

Thin, long
elongation tip
for easy root
canal application

Game-changing syringe.

Proprietary 3M design takes automix cement applications to the next level. Developed with your needs in mind for easy dispensing, cleaning and storage. The small, low-waste 3.4g syringe provides enough material for between 15–20 applications.







Easy excess clean-up.

Efficient to the last step. Its unique chemistry ensures excess outflow stays put at the restoration margins, and can be easily removed after a 2–3 second tack cure.



Self-adhesive cementation of a 3M™ Lava™ Esthetic Fluorescent Full-Contour Zirconia bridge; Image courtesy of Dr. Jan Gueth, Germany

So simple yet so powerful.

-  **Two components for virtually all indirect and direct bonding indications**
-  **Outstanding self-adhesion of RelyX Universal Resin Cement**
-  **Scotchbond Universal Plus Adhesive further enhances bond strength to all substrates**
-  **Unique 3M self-sealing syringe**
-  **Easy excess clean-up**
-  **High aesthetics**




The strong partner for adhesive cases.

3M™ Scotchbond™ Universal Plus Adhesive is the first radiopaque, all-in-one universal adhesive for all direct and indirect bonding indications. And, it is a fully aligned system: the adhesive is cured by the cement, so no light-cure is needed.


Scotchbond Universal Plus Adhesive serves as universal primer for all materials including glass ceramics and further boosts the cement's excellent self-adhesive bond e.g. to zirconia.

Outstanding self-adhesion for maximised efficiency.

With its superior self-adhesive bond strength to dentin, RelyX Universal Resin Cement reliably treats most cases with no need for an additional adhesive or primer. For light- or self-cure.

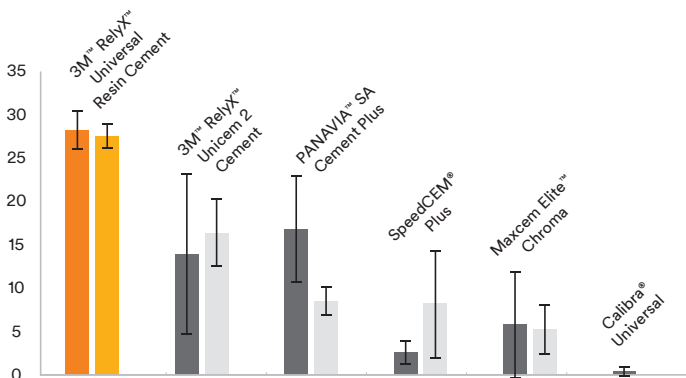

Shear bond strength (MPa) to dentin after artificial aging
R. Afutu, M. Abreu, G. Kugel, Tufts University, J. Dent. Res. Vol 98A, No 3629, 2019

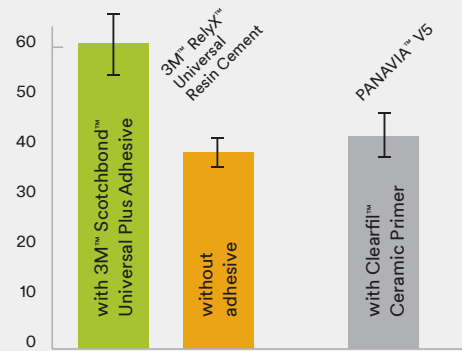

Less stock, more savings

95%

of dentists agreed the universal character of RelyX Universal Resin Cement helps to **reduce stock and save costs.**

Easier for operator staff

97%


of dentists agreed using one universal versus multiple resin cements is **easier for operator staff.**

Virtually no post-operative sensitivity

99%

of dentists reported that RelyX Universal Resin Cement shows **virtually no post-operative sensitivity*****

High aesthetics with four colour-stable, fluorescent shades

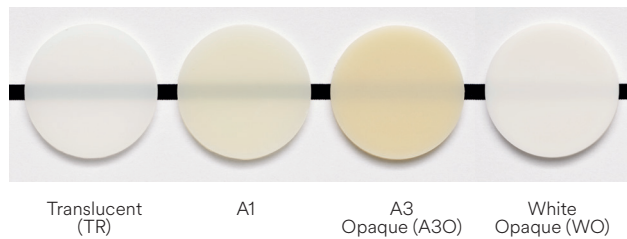

Ordering information

3M™ RelyX™ Universal Resin Cement

Item no. Product information

Intro Kit

56968 4x 3M™ RelyX™ Universal Resin Cement syringes shades TR, A1, A3O, WO (3.4g each), 60x 3M™ RelyX™ Universal Micro Mixing Tips, 20 Elongation Tips, 1 vial 3M™ Scotchbond™ Universal Plus Adhesive (5 ml), 1x 3M™ Scotchbond™ Universal Etchant (3 ml), 25 blue etchant tips, 50 orange microbrushes, step-by-step guides

Trial Kits

Each trial kit contains: 1x 3M™ RelyX™ Universal Resin Cement syringe (3.4g), 15x 3M™ RelyX™ Universal Micro Mixing Tips, 5 Elongation Tips, 1 vial 3M™ Scotchbond™ Universal Plus Adhesive (1.5 ml), 1x 3M™ Scotchbond™ Universal Etchant (3 ml), 25 blue etchant tips, 50 orange microbrushes, step-by-step guides

56969 Shade TR 56970 Shade A1

Refills

Each refill contains: 1x 3M™ RelyX™ Universal Resin Cement syringe (3.4g), 15x 3M™ RelyX™ Universal Micro Mixing Tips, 5 Elongation Tips

56971 Shade TR 56973 Shade A3O

56972 Shade A1 56974 Shade WO

Value Packs

Each value pack contains: 3x 3M™ RelyX™ Universal Resin Cement syringes (3.4g each), 45x 3M™ RelyX™ Universal Micro Mixing Tips, 15 Elongation Tips

56977 Shade TR 56978 Shade A1

Mixing Tips

56975 30x 3M™ RelyX™ Universal Micro Mixing Tips

56976 15x 3M™ RelyX™ Universal Micro Mixing Tips with Elongation Tips

3M™ Scotchbond™ Universal Plus Adhesive

Vial

41293 **Intro Kit:** 1x 3M™ Scotchbond™ Universal Plus Adhesive vial (5 ml), 50 disposable brushes, 1x 3M™ Scotchbond™ Universal Etchant (3ml), 25 dispensing tips

Refill: 3M™ Scotchbond™ Universal Plus Adhesive vial (5 ml)

41294 1x 5 ml vial 41295 3x 5 ml vial

Unit Dose

3M™ Scotchbond™ Universal Plus Adhesive Unit Dose L-Pops, 1x 3M™ Scotchbond™ Universal Etchant (3ml), 25 dispensing tips

41297 **Intro Kit: 50x L-Pop 41299 200x L-Pop**

41298 100x L-Pop 41304 400x L-Pop

3M™ Scotchbond™ Universal Etchant Etching Gel

41263 2x 3M™ Scotchbond™ Universal Etchant (3ml each), 50 dispensing tips



3M™ RelyX™ Universal Resin Cement Intro Kit (56968)



3M™ RelyX™ Universal Resin Cement Trial Kit TR (56969)



3M™ RelyX™ Universal Resin Cement Refill TR (56971)



3M™ Scotchbond™ Universal Plus Adhesive Refill (41294)



3M™ Scotchbond™ Universal Plus Adhesive Unit Dose 100 L-Pop (41298)

3M™ Scotchbond™ Universal Plus Adhesive Unit Dose:
For efficient hygiene management. Reduces the risk of cross-contamination.



3M Oral Care UK
47 Charnwood Campus
Bakewell Road
Loughborough
LE11 5RB

3M.co.uk/dental

(+44) 0870 536 036

3MOralCareUK@mmm.com

@3MOralCareUK

3M, Lava, RelyX and Scotchbond are trademarks of 3M Company or 3M Deutschland GmbH. Used under license in Canada. All other trademarks are owned by other companies.

© 3M 2020. All rights reserved. OMG119934.